



HPV/p16 Analyte Control

Cell line controls for immunohistochemistry
and in situ hybridization.

Research Use Only

PRODUCT AVAILABILITY

Product Code	Product Description
HCL004	X2 Cut slides
HCL005	X5 Cut slides
HCL006	X1 Cell microarray block

APPLICATION

This product is suitable for use in immunohistochemistry and in situ hybridization.

MATERIALS

Three formalin fixed paraffin embedded cell lines. Two of which contain incorporated HPV 16 and 18 viral genes.


Cell line A: Human breast adenocarcinoma

Cell line B: Human cervical adenocarcinoma


Cell line C: Human epidermoid carcinoma

Cells are fixed in 10% neutral buffered formalin and paraffin wax embedded.

Sections are cut at 4µm, mounted on positively charged slides and baked overnight at 37°C.

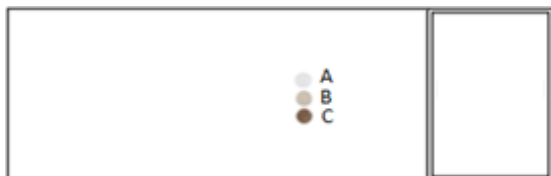

Cell microarrays (CMA) contain cores that are 1.5-2mm in diameter and 3-3.5mm in length. It is possible to obtain over 300 sections depending on thickness.

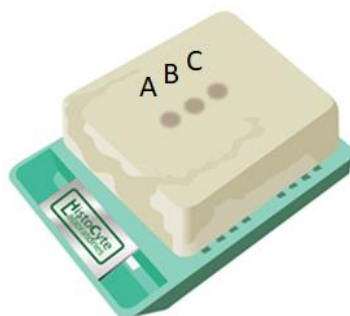

Expression Profile

Cell Line	HPV Gene Copy	mRNA E6/E7 Copy	p16 Expression
A	Negative	Negative	Negative
B	Medium	High	High
C	High	High	High

Storage and Handling

Store at 2-8°C. Do not freeze (for expiration date please see the product label)

WARNINGS AND PRECAUTIONS

1. The product is intended for research use only.
2. It is the responsibility of the end user to determine suitability with their reagents and procedures within their laboratory.
3. Do not use after expiration date printed on product labels. The user must validate any storage conditions other than those specified in the package insert.

TROUBLE SHOOTING

For further help please feel free to contact HistoCyte Laboratories Ltd at info@histocyte.com or call +44 (0)191 603 1007.



For updates and additional product information please visit: www.HistoCyte.com